BIOS living

# Living<sup>2</sup> Lift Chair

---

## LF900

Use and care manual

# Dear Customer

Thank you for purchasing LF900 recliner chair. Please read this book carefully before using armchair.

## SUMMARY

Safety Precautions.....	03
Product Installation.....	05
Moving Product.....	07
Grounding Instructions.....	07
Product Characteristics.....	08
Remote Control Use.....	09
Troubleshooting.....	10

# Safety Precautions

Some precautions must be taken when using LF900 recliner chair in order to avoid any risk of injury. Please read carefully the followings instructions carefully.

 **DANGER:** To reduce the risk of electrical shock:

- Always unplug this product immediately after using and before cleaning.

 **WARNING:** To reduce the risk of burns, fire, electric, shock or injury to persons:

- Do not put hands or fingers near the mechanisms
- If the power cable is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard.
- Ensure that the chair is connected to a proper electric installation.
- Always unplug the chair after use and before cleaning.
- Do not use the chair in a humid room.
- Do not expose the chair to direct sun light.
- In case of blackout, please turn power switch OFF in order to avoid any damage to the product when power returns.
- Never plug or unplug the chair with wet hands.
- Close supervision is required when the chair is used by children, invalids or disabled persons.
- No more than one person should use the recliner chair at the same time.
- Use the chair only for the functions explained in this manual and do not use accessories not supplied by the manufacturer.
- Do not try to stand up on chair's seat or arm because you could damage product or have an accident.
- Do not sit on the remote control.
- Avoid falling asleep while using the recliner chair.
- Always keep the chair cleaned and never insert objects into gaps.
- Do not move the chair by pulling power cord.
- Do not pull on the AC lead.
- Ensure that no one, no object or pet is under the chair before reclining it.
- This product is designed for household use only.
- Make sure that upholstery is in good condition before using the chair.
- When using this product after it as been in storage, check to see that the product is operating properly prior to use.
- Do not use the chair in a room where aerosol are sprayed and where oxygen is dispensed.

## SAVE THESE INSTRUCTIONS

**Do not attempt to repair this product yourself.**

This product is designed for home use. Commercial use will void warranty.

### Correct Disposal of this product

This marking indicates that this product should not be disposed of with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmentally safe recycling.


**Note:** Our firm has the right to re-design and improve the product. Should you encounter some differences, we apologize as the product could be different from that in your manual.

# Installation


---

**Note :** Do not use sharp or pointed objects to remove chair's packaging. Such tools can damage the upholstery of the product.

**Step 1 :** Make sure there is at least 13cm of free space behind the backrest in raised position and 63cm in front of the leg rest in folded position when arranging the armchair.

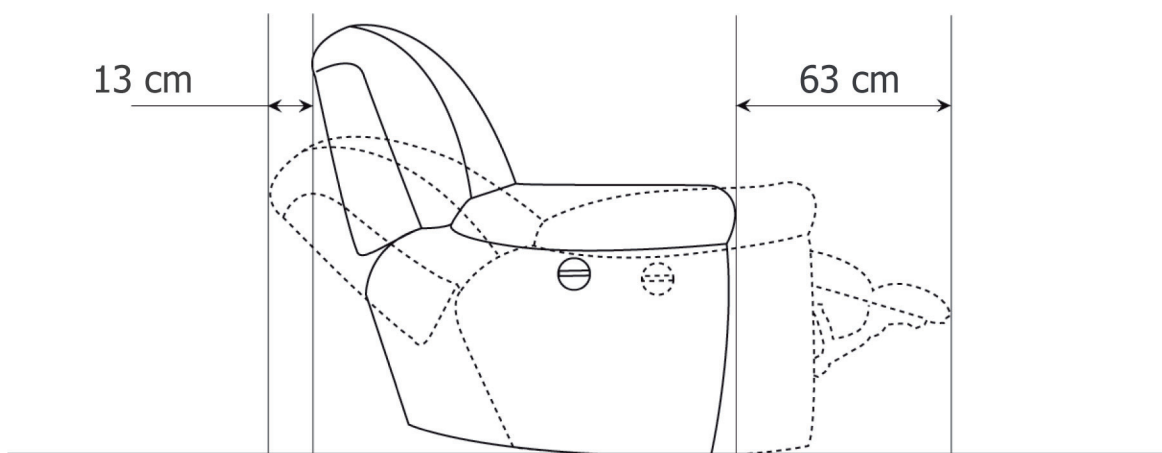
**Step 2 :** Make sure the backrest reclining system is in raised position, then insert left and right flat head bolts of seat (d) into left and right backrest's brackets (c).

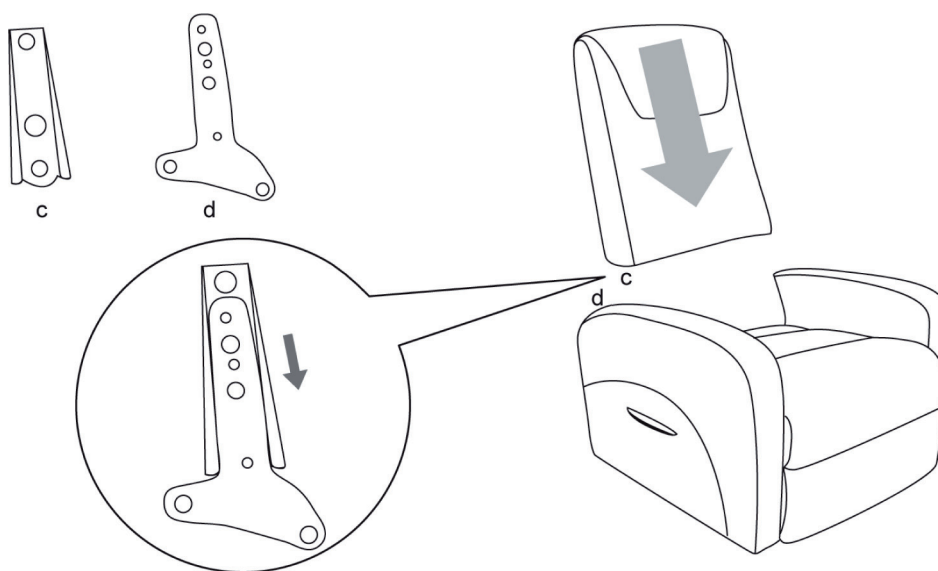
### Step 3 : Power cord connection

- Connect the three electric cords as shown in the diagram below (A & B)


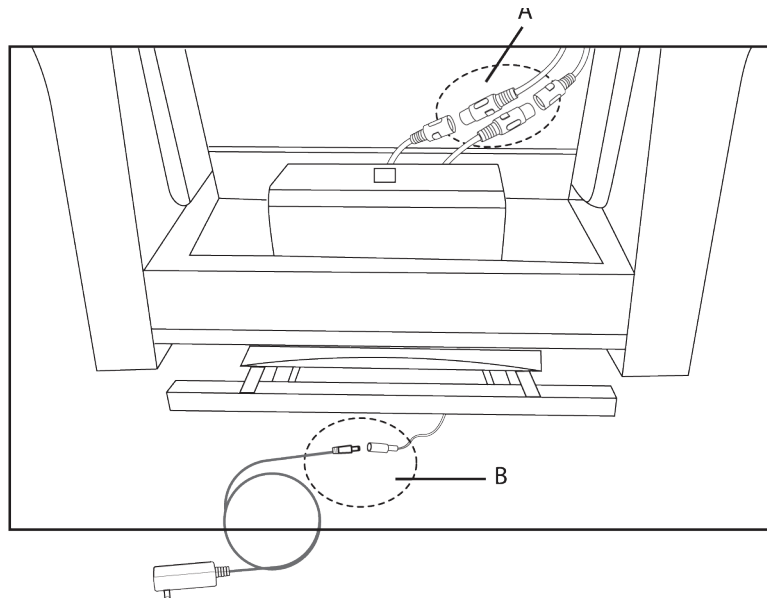
- It is strongly suggest you put two 9 V (DC) batteries into the compartment located under armchair's seat. This allows user to use reclining system in case of a blackout. (Batteries NOT included)
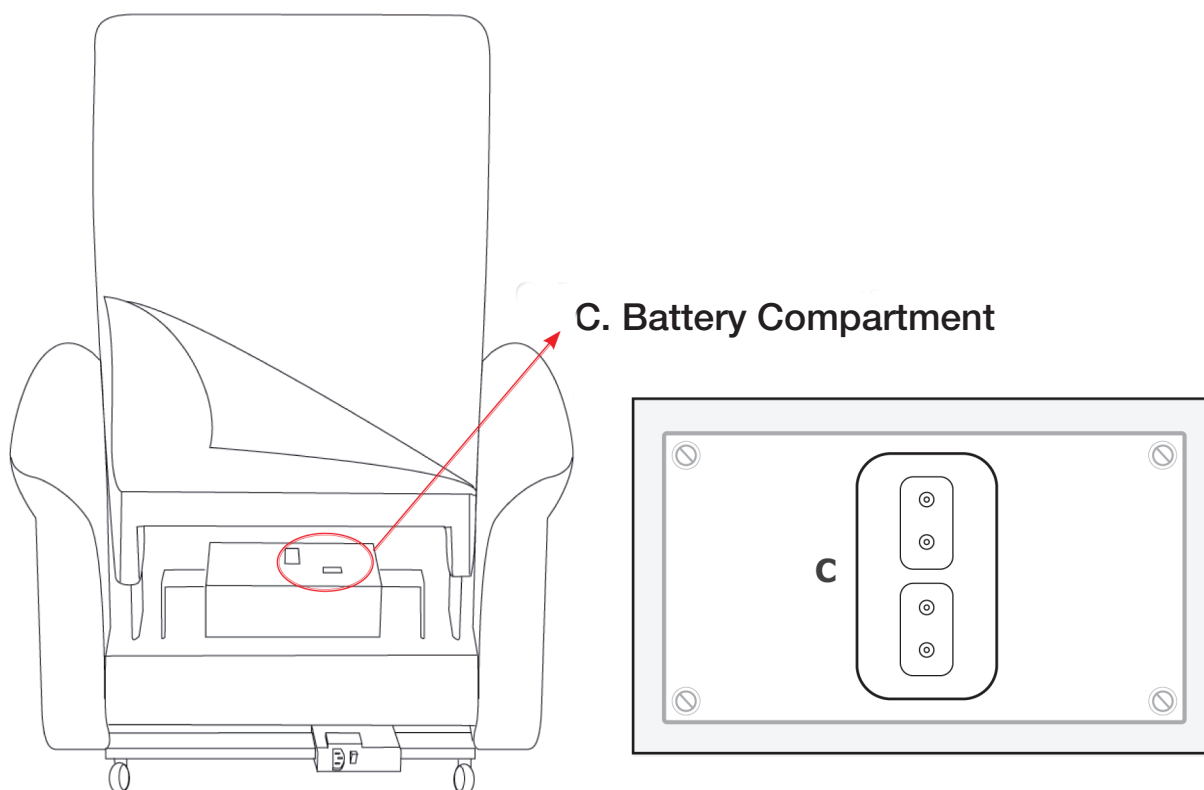
# Moving

Raise up the front part of the chair by holding side panels and use the wheels located on rear part of chair's base in order to move recliner chair.

**Note:** Do not hold armrests to move the chair. Do not drag or push the chair on fragile surfaces.

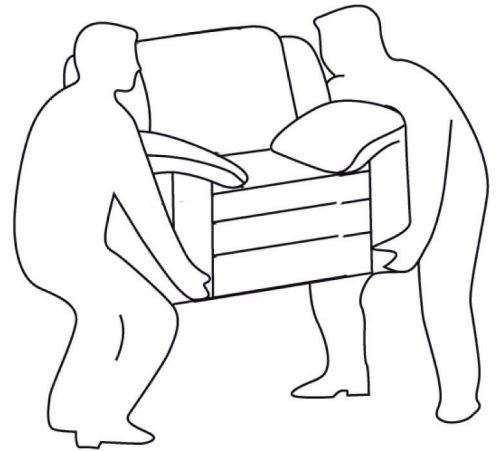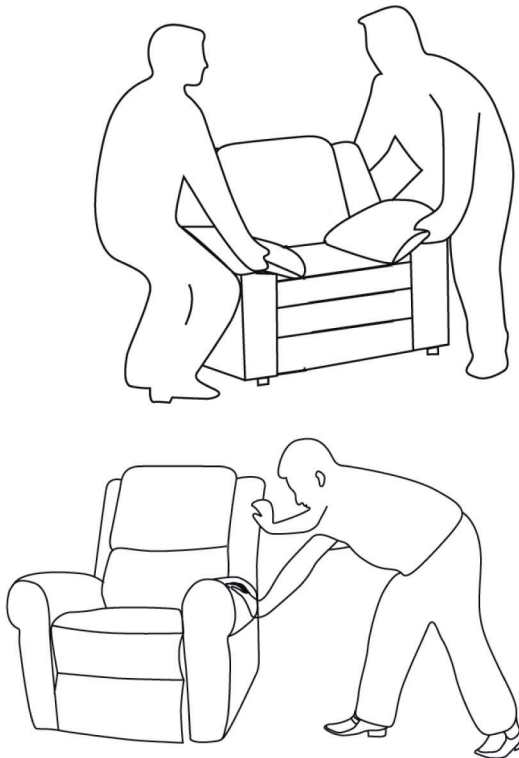
## Grounding instructions

Insert power plug into an appropriate outlet.

Make sure that electrical installation is regular in order to avoid any electrical shock. If not, please contact a qualified electrician.

# Product Characteristics


Information also available on the information tag on armchair's main frame.
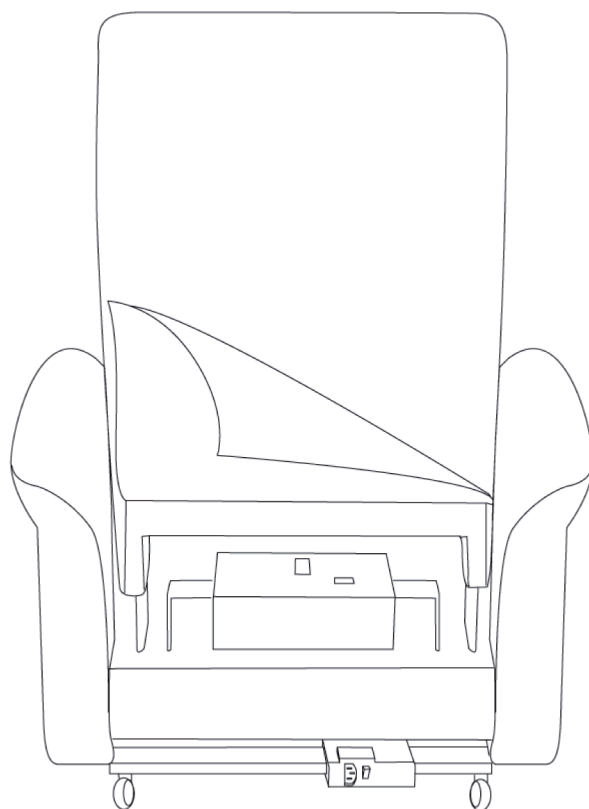
<b>Model No :</b>	<b>ES-8231</b>
<b>Weight :</b>	<b>54 kg</b>
<b>Voltage :</b>	<b>120 V</b>
<b>Frequency :</b>	<b>60 Hz</b>
<b>Input :</b>	<b>24V 3A</b>
<b>Power Consumption :</b>	<b>72 W</b>
<b>Certification :</b>	<b>UL</b>

# Remote Control Use
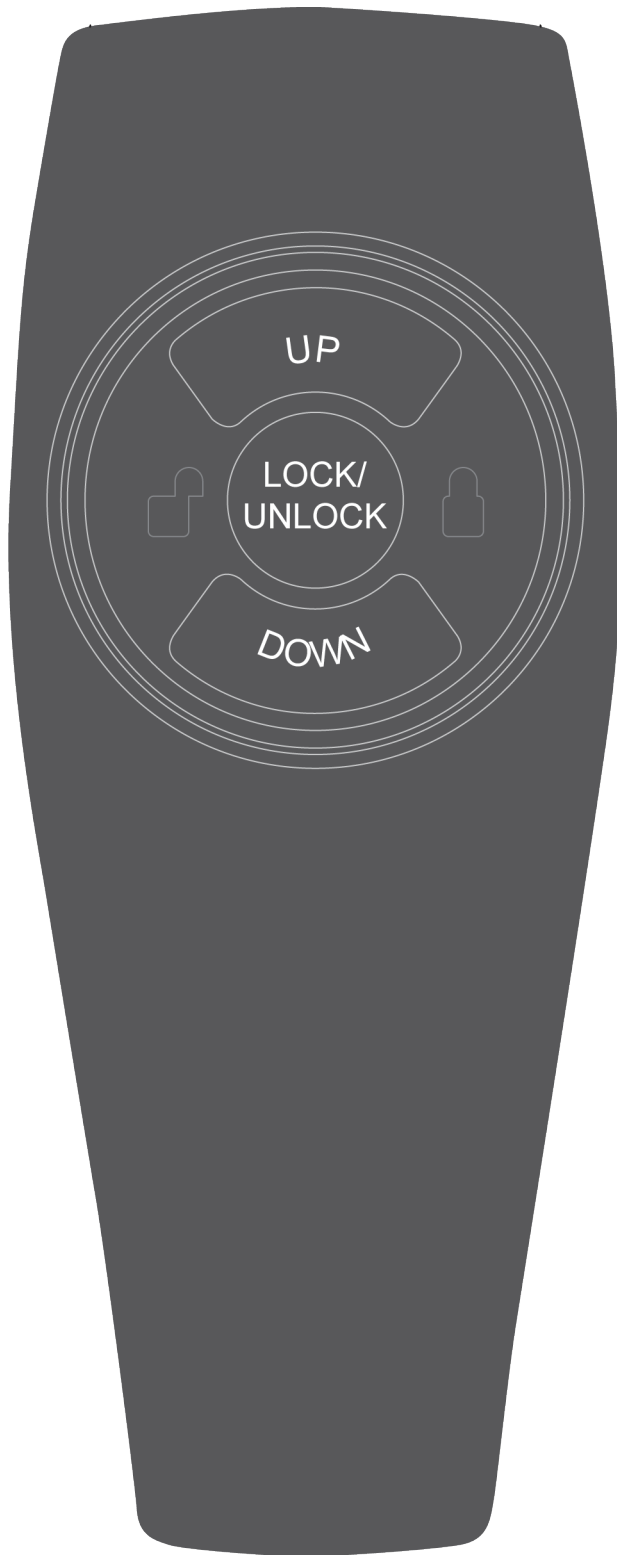
## 1. UP

Press and hold this button to raise up backrest and recline leg rest. When backrest is already in vertical position, use this button to lift up chair's seat. Release this button when you have reached to desired position.

## 2. LOCK/UNLOCK

Press this button to lock your position selection in place. Press it again to unlock armchair and regain control.

## 3. DOWN

Press and hold this button to raise up leg rest and recline backrest. When seat is in upper position, then use this button to move seat down.

# Troubleshooting

1. Recliner chair does not work even when power is ON.

**Actions:** Check that power plug is properly connected to power socket.

2. Remote control is ON but chair does not to work.

**Cause:** After a prolonged use, chair can be overheated.

**Solutions:** Please switch OFF the chair and let it cool for 30-50 minutes.

Please contact your retailer if problems persist.

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